

## Vertical Autoclave SLES Series


Cominus Scientific makes for reasearch work. We are the most trustworthy and reputable manufacturers of cautoclave. We work hard to deliver high-quality responses with comparable performance par excellence for our customers at the lowest possible price.

With the goal of boosting laboratory productivity, the Cominus vertical floor-standing laboratory autoclaves with top-loading access meet the essential requirements for general labware sterilients, and research institutions in numerous enterprises, educational and medical establishment

They are designed to sterilise a variety of solid and liquid materials, including waste bags, glassware, plastic, metallic equipment, culture medium, and other laboratory supplies.

### Special Features:-

- **ECONOMIC AND DURABLE:-**
  - SLES Series autoclaves are economic and robust autoclaves with excellent performance for general laboratory sterilization procedures. They can be used either for solids and liquids sterilization procedures and they consume limited valuable laboratory resources such as water, power or operator time.
- **EASY INSTALLATION AND MAINTENANCE**
  - Every SLES Series autoclave is a plug and play equipment that do not need dedicated installation connections. They simply need a power source and can work even without a connection to the drainage. All models include casters so they can be used in different areas of the same facility
- **SAFETY FIRST**
  - SLES Series autoclaves are equipped with several features to ensure the safety of the operators. These include an overpressure safety valve, a thermally insulated door, an overtemperature safety thermostat, an open door detection system and an independent safety pneumatic system that locks the main door while positive pressure exists inside the sterilization chamber.

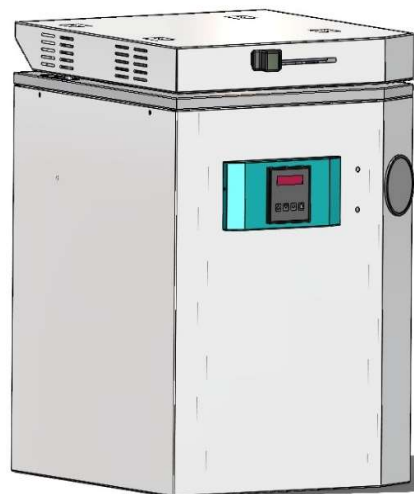



### **Technical Specifications**

<b>Model No</b>		<i>CM-VAV-SLES-50</i>
Control System		PID Microcontroller
Display		LCD
Internal Chamber Volume	Ltrs	50
Chamber Type		Vertical
Door Opening Mechanism		Sliding
Operating Temperature	°C	100-125
Operating Pressure	PSIG	15
Castor Wheel		Yes
Operation		Fully Automatic
Power		220 V AC 50 Hz
Warranty		1 Year

*Check for the latest cominus products*

<https://cominus.in/>

**Office:- D-132 Sec 4 Bawana  
DSIIDC Complex Delhi-110039**

**Landline:- 011-47553098**

**Mob:- 9027822747, 9667529087**

**Email:- sales@cominus.in**